

Contour and straight shape

Silicone bladder, easy to wipe - easy to dry

Steam autoclavable at 134°C

UDI code and serial number



## Tourniquet Silicone Cuff

Tourniquet cuff, reusable through automatic processing

# Tourniquet Silicone Cuff

Tourniquet Silicone Cuff is a reusable tourniquet cuff made of durable, heavy-duty material and is suitable for repeated use after sterilization by steam autoclave at a temperature of 134°C for 5 minutes. The cuff is available in straight shape (for small or midsize upper extremity) and contour shape (mainly for lower conical shaped extremity) and is provided in different sizes, which are colour coded for optimal cuff selection.

## Features




### Cuff Design

- Large hook and loop fastener for reliable cuff fixation
- Colour coding makes size selection and ordering easier
- Tubing made of kink-resistant silicone material
- Positive Locking Connector made of metal including Sealing Cap
- Contour and straight shape
- PVC-free




### Special Features

- Silicone bladder, easy to wipe - easy to dry
- UDI code for product traceability
- Serial number for tracking of product usage
- Processable up to 50 times (steam autoclave)

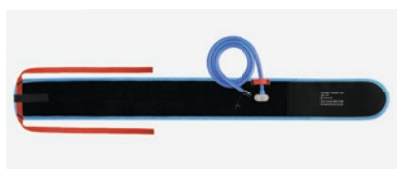


### Straight Shape

- Cylindrical design
- Does not show any difference in diameter between proximal and distal part

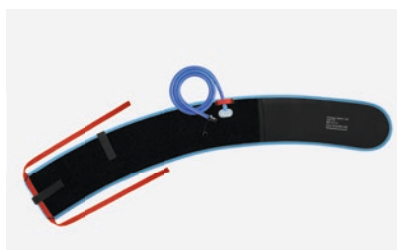

### Contour Shape

- Conical design with a smaller diameter distally than proximally when applied
- Occludes the blood flow at lower cuff pressures due to an improved fit to the extremity
- Specific bladder design (wave pattern) to prevent cuff migration

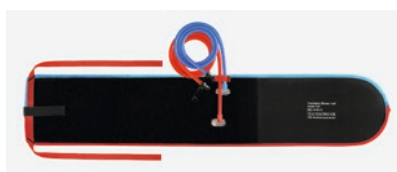

### Double Cuff for Intravenous Regional Anaesthesia (IVRA)

- Double cuff with 2 individually inflatable bladders

## Order information

### Tourniquet Silicone Cuff / Single cuff, reusable, autoclavable

Length	Width	Colour	REF	Box
20 cm	2.5 cm	-	20-64-700	1
30 cm	4.5 cm	■	20-64-710	1
35 cm	7.0 cm	■	20-64-611	1
46 cm	8.5 cm	■	20-64-612	1
46 cm	11.0 cm	■	20-64-512*	1
61 cm	11.0 cm	■	20-64-522*	1
76 cm	11.0 cm	■	20-64-527*	1
86 cm	11.0 cm	■	20-64-528*	1

### Tourniquet Silicone Cuff / Double cuff, reusable, autoclavable

Length	Width	Colour	REF	Box
30 cm	8 cm	■	20-60-710	1
46 cm	12.5 cm	■	20-60-712	1
61 cm	12.5 cm	■	20-60-722	1

\*Contour shape

### Sealing Cap / For Tourniquet Silicone Cuff, reusable

REF	Box
22-60-111	5

This product is manufactured without the use of natural rubber latex, unless otherwise specified.  
This product does not contain phthalates which require labelling according to CLP Regulation (EC) 1272/2008.